













# **GENERAL SPECIFICATIONS** (I.on/I.on PRO)

- Dimensions:
- 225 mm (W)
- 225 mm (H). - 69 mm (D).
- · Weight:
- 1.2 Kg (basic) to 1,9 Kg (complete).
- Non-rechargeable internal battery:
  Type: Manganese Lithium Dioxide (LiMnO<sub>2</sub>) 18V,
- 1400 mAh (optional 1) or 2800 mAh (optional 2).
- Duration of optional battery 1 (1400 mAh): more than 80 shocks in 200 J or 6 hours of continuous monitoring.
- Duration of optional battery 2 (2800 mAh): more than 300 shocks in 200 J or 15 hours of continuous monitoring

# Rechargeable internal battery:

- Type: Li-lon, 14.4 VDC 4.0 A/h.
- Duration: 18 hours of continuous monitoring.
- Time for battery full charge (fully discharged): 5

# • Battery power supply charger: - Power 100 - 220 V/50-60 Hz.

- Consumption (maximum): mains supply 1 A.
- Output: 16.8 VDC. 1 A.
- Use only source manufactured by Instramed.

#### · Battery storage:

- Storing the battery for a long period of time in temperatures higher than 35°C (95°F) will reduce its capacity and shelf life.

# • Pre-adjusted defibrillation scales:

- Adult (automatic): 1° shock 150 J; following shocks 200 J.
- Child: 50 J.
- · Defibrillation scales adjusted by user (via SoftDEA):
- Adult (non-rechargeable battery 1400 mAh): scales between 120 J and 200 J.
- Adult (non rechargeable battery 2800 mAh or rechargeable): scales between 120 J and 360 J.

- Internal memory storage:
  100 events or 2 hours of ECG recording.
- · Ambient sound storage:
- Up to 10 hours (optional)
- Protection rating:
- · Classification:
- Class II, internally energized equipment.

#### · Electrical insulation: CF type

- · Operating mode: Continuous operation
- · Maximum time from rhythm analysis beginning to discharge readiness:
- 200 J: 20 seconds.
- 360 J: 25 seconds
- · Maximum time from start of operation beginning to full discharge readiness: - 200 J: 25 seconds.
- 360 J: 35 seconds.

#### Non-frequent use equipment:

- Comply with the requirements for non-frequent use equipment as specified in NBR IEC 60601-2-4

#### **ENVIRONMENTAL SPECIFICATIONS**

- Temperature:
- Operational: 0 to 50°C. Storage: -20 to 50°C.

# • Humidity:

- Operational: 10 to 95% RH, without condensation.
- Storage: 10 to 95% RH, without condensation.

# · Altitude:

- Recommended to use in altitudes lower than 2000 meters.

## **DEFIBRILLATOR**

#### · Waveform:

- Biphasic truncated exponential. Wave shaped parameters adjusted according to the patient's impedance.
- · Shock application:
- By means of multifunctional adhesive pads.
- · Commands:
- Front panel button: on/off.

# Scales for defibrillation:

- Adult: adjustable from 120 up to 360 J (depending on the battery and via SoftDEA). 1° shock - 150 J; following shocks - 200 J.
- Child: 50 J.
- · Adults/children selection:
- Automatic due to the size of the pads.
- · Charge command:
- Automatic after identifying arrhythmias that should receive shock.
- Shock command:
- Front panel button, when blinking.

# · Maximum charging time:

# Rechargeable battery

- 50 J: < 2 seconds
- 150 J: < 3 seconds. 200 J: < 4 seconds.</li>
- 270 J: < 5 seconds.</li>
- 360 J: < 6 seconds.

## Non-rechargeable battery (1400 mAh):

- 50 J: < 3 seconds.
- 150 J: < 8 seconds.
- 200 J: < 11 seconds.

#### Non-rechargeable battery (2800 mAh):

- 50 J: < 2 seconds.</li>
- 150 J: < 5 seconds.
- 200 J: < 6 seconds.
- 270 J: < 8 seconds. - 360 J: < 10 seconds.

# · Pads size:

- Adult = area: 82 cm<sup>2</sup> (32.3 in<sup>2</sup>).
- Child = area: 30 cm<sup>2</sup> (11.8 in<sup>2</sup>).
- · Pads cable length: - 2 meters.

# Maximum output voltage:

- · Maximum output current: - 80 A (25 Ω).

#### OTHER SPECIFICATIONS (I.on PRO)

#### Commands:

Touchscreen: allows for the selection of manual mode and defining energy scales.

# · Scales for defibrillation (Manual mode)

- Adult: scales between 120 J to 360 J.
- Child: 10, 20, 30, 40 or 50 J.

# **STANDARDS**

- IEC 60601-1
- IEC 60601-1-2
- IEC 60601-1-6 - IEC 60601-1-8
- IEC 60601-1-9
- IEC 60601-2-4 or equivalent IECs.

\*Some items are optional.

# Instramed Indústria Médico-Hospitalar Ltda.

Beco José Paris, 339 - Pavilhão 19 Tel.: +55 (51) 3073 8200

CEP 91140-310 Porto Alegre RS Brasil


E-mail: comex@instramed.com.br WWW.INSTRAMED.COM.BR

European Representative:

## Obelis S.A.

Bd. Général Wahis 53, 1030 - Brussels, Belgium Phone: + 32.2.732.59.54 | Fax: + 32.2.732.60.03 E-mail: mail@obelis.net



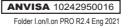



# INSTRAMED





\*Some items are optional.

- Much more economical disposable pads (and completely independent of battery replacement and/or the CPR feedback accessory).
- SoftDEA software included for connection, feature settings, data downloading and viewing using a PC (via USB).
- · Easily switch between automatic and manual modes with a simple touch on the screen (I.on

SIMPLE TO USE

• It evaluates, with

variables and

shock therapy.

sophisticated sensors,

the patient's condition.

considers the clinical

automatically applies

the most appropriate

• It allows any individual

with basic training to

perform the care of

a victim in cardiac

arrest, facilitating

and multiplying the

possibilities of life-

• Just press the single

front button and follow

the guidance by voice

and visual indicators.

saving.

**NEW!** MANUAL **SELECTION BUTTON** ADULT/ CHILD

> It allows the use of regular pads, even in child patients.

> > **ADVANCED**

• The I.on PRO model

manual operation,

thus allowing health

professionals, based

and in the ECG curve

screen, to personalize

on their experience

shown in the device

the parameters to

be applied in the

• In the touch screen

with excellent contrast

and viewing area, the

operation mode and

then the charge up to

user selects the manual

treatment.

360 J.


adds the flexibility of




Positioned on the patient's chest, the CPR Maestro (optional) provides real-time feedback on the rescuer's performance by guiding him through voice and LCD display about the procedure's optimal frequency, strength and interval

The result is a much more effective and consistent performance of CPR, resulting not only in increased chances of patient survival, but reducing the occurrence of sequelae due to lack of oxygenation



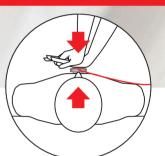



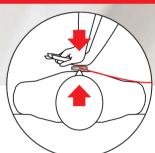



 Via the SoftDEA application (included), the Auto-Charge Sequence function can charge preset energies for the first, second and third shocks.

INTELLIGENT

- Through the 3-lead ECG cable, the I.on can act as an **Automated External** Defibrillator or ECG monitor, alarming when it detects a cardiac arrest situation.
- It has built-in microphone and internal memory allowing digital storage of up to 10 hours of ambient sound (optional).

- Instramed technology with robust design used in hundreds of medical organizations, developed for the reality of emergency care.
- · A wide network of distributors providing sales of products and parts, along with qualified technical assistance, that provides excellent postsales service.
- With more than 30 vears in the business. Instramed is a manufacturer that maintains stock parts and components even for discontinued products.